




A Simple and Effective Training Unit

# iFMT - 2001D

INTELLIGENT FOOT MUSCLE TRAINER

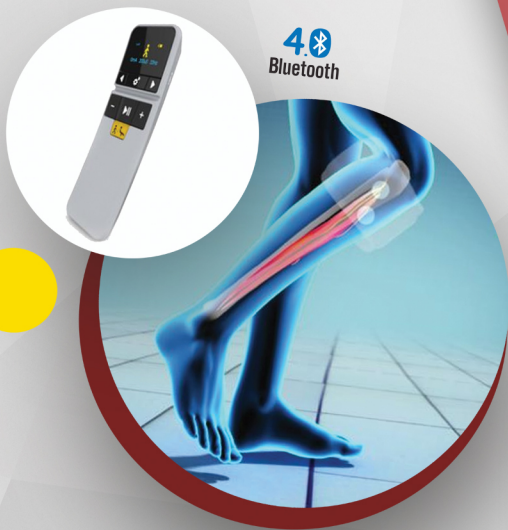

iFMT - 2001D can be used as a strengthening exercising and walking aid for users with weak foot and leg muscles.

## INNOVATIONS :

- Advanced ergonomic design makes the electrode fit the leg muscles.
- Ultra-thin design makes the device less noticeable and can be covered with trousers perfectly.
- Magnetic clasp for one-handed operation.

## MODES

TRAINING MODE

GAIT MODE

## WORKING PRINCIPLE :

iFMD-2001D Intelligent Foot Muscle Trainer adopts advanced sensor technology and intelligent algorithms, precisely controlling the time and duration of stimulation by tracking the swing angle and pace of user's legs. Those mild pulses stimulates patients' leg muscles, making them lift the affected foot at an appropriate phase, while walking and therefore enabling patients to walk more steadily, naturally and safely.